

King Single Mattress Replacement System




FEATURES

- ★ Stage II medical mattress replacement system for medium pressure care treatment
- ★ A-B lines air tubes alternating for better body / pressure relief
- ★ Top grade TPU coated with Nylon material
- ★ 4-way stretch PU cover with zipper closure, fire retardant and waterproof
- ★ Pillow function in static to provide redundancy through power failure or air cell damage
- ★ CPR release for quick air release in emergency
- ★ Warranty: 1 year

SPECIFICATIONS

- ★ Size: 80"L*42"W*9.5"H King Single Size Mattress Fit
- ★ 4.5"H 18 cells and 5"H 16 cells for more options(Nylon TPU)
- ★ Mattress material: Nylon TPU
- ★ Cover material: PU
- ★ Bottom material: Dacron PVC
- ★ Net weight: 8KG
- ★ Max Weight Capacity: 145KG

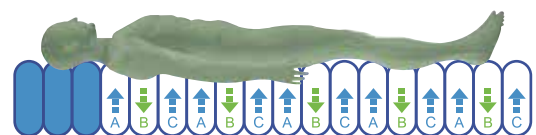

Air Pump Unit

FEATURES

- ★ Aluminum motor cover, excellent heat dissipation, anti-vibration
- ★ Ultra silence, noise only ≤20dB
- ★ Digital LED display screen showing weight and pressure
- ★ Auto pressure control, static mode, sound alarm
- ★ Lock mode, seating mode
- ★ Warranty: 1 years

SPECIFICATIONS

- ★ Power: 220V-240V / 50Hz
- ★ Air output: ≥6L/min
- ★ Pressure range: 14-62mmHg
- ★ Cycle time:
 - 14mins(A-B lines alternation)
 - 21mins(A-B-C lines alternation)
- ★ Power consumption: 15VA
- ★ Fuse: 0.5A
- ★ Size: 275*135*94mm(L*W*H)
- ★ Net Weight: 1.3KG

